106x150mm

22.4.8

DIGITALES THERMOMETER KFT-04

! WARNUNG

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Fabricante/Manufacturer:
Kangfu Medical Equipment Factory
N°380 Ningkang East Road, 325600 Yueqing,
Zhejiang, China

EC REP Luxus Lebenswelt GmbH
Kochstr. 1, 47877, Willich, Germany

CE 0123

106x150mm

22.4.8

User manual in Deutsch,
Français, Portugais, Italiano.

2. Boost up slightly, then clamp it, in order to full touch between thermometer and body, please hold the arm lightly

Boost up slightly
Clamp
Hold the arm Lightly by hand

DIGITAL THERMOMETER

KFT-04

WARNING

- Don't be auto diagnosis and treatment according to result of body temperature, consult your doctor, auto diagnosis may cause the worsening of your illness.
- This product is used to your **oxter**, mouth and rectum. Please don't use to other places (E.ear ETC.) Otherwise unable to get correct body temperature, even hurt your ear ETC.
- Please place it to untouched for children, don't use it alone by children, otherwise it will be hurted to your children after use.

ATTENTION


- This product (IPX0 device) is no waterproof, have to avoid water ETC. Liquid to come into the product, otherwise the measured value may cause error or breakdown.
- Please don't bite this product, otherwise cause to occurrence and damage.
- Please put battery, screw to untouched place for children, please go to doctor if children accidentally swallowed it.
- Please wash with plenty of water in case the liquid from battery into your eyes, you must go to doctor for treatment, otherwise you may loss your sight ETC. **disability**.
- Please don't throw the battery into the fire, otherwise the battery will be broken.
- Please don't use it at the place of strong static electricity or electromagnetic wave, otherwise the measured value may cause error or breakdown.

METHOD OF CORRECT MEASUREMENT

The measured value may be get incorrect of not proper operation when you use it.

WHEN MEASURE AT OXTER Highest temperature at the middle of oxter

- Attach to the middle of oxter

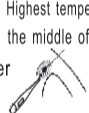WRONG MEASUREMENT METHOD AS FOLLOWS

Thermometer is not in proper place Thermal is easy to be out

WHEN MEASURE IN MOUTH

- Attach to the right/left of tongue's root.

- To coarctate the thermal of thermometer by tongue

- Carry it by hand to avoid the movement of thermometer, don't open your mouth during measurement.

WHEN MEASURE AT RECTUM

This is most accurate way to measure, and specially for baby and KID. Please insert the thermal 2-3cm deep to rectum.

THE FOLLOWING SITUATIONS NOT BE GOT CORRECT TEMPERATURE

- Just finished the sports, bathing or eating;
- In the quilt very long time, sweat in the oxter;
- Exercise after getting up;

USAGE METHOD

- Press "power" button, after the beep, the screen is shown as right FIG. **188.8°C** or **188.8°F**
- Thermometer display the temperature last time. **37.0°C** or **98.6°F**


- When it is ready to measure, put thermometer in correct place after "°C" or "°F". 
- When basically the temperature is no variation, °C or °F stopped flashing, and buzzer makes a sound, remind you to confirm the result to measurement.
- Please keep the thermometer as long as possible to get accurate temperature
- Press "power" button again, switch off.
- Thermometer is self-closing after 10 min. In order to get accurate temperature 2min at oxter 80 seconds in the mouth 60 seconds in rectum
- Press and hold "power" button 3 seconds to switch between °C and °F

METHOD BATTERY REPLACEMENT

If the screen is shown as right **37.0°C** or **98.6°F** FIG. Please replace the battery.

Method

- Open battery cover
- Push battery out by small bar as per arrow shown below
- Push battery in the direction of arrow, -positive face up

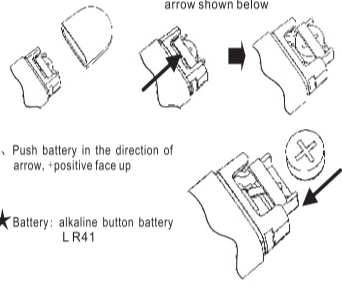★ Battery: alkaline button battery LR41

Flashing

DISPOSAL PROCEDURE(DIR.2002/96/EC-WEEE)

The symbol on the bottom of the device indicates the separated collection of electric and electronic equipment. At the end of life of the device, do not dispose it as mixed solid municipal waste, but dispose it referring to a specific collection centre located in Your area or returning it to the distributor, when buying a new device of the same type to be used with the same functions. devices is carried out forecasting a European environmental policy aiming at safeguarding, protecting and improving environment quality, as well as avoiding potential effects on human health due to the presence of hazardous substances in such equipment or to an improper use of the same or of parts of the same.

WARRANTY

This appliance is warranted for a period of 1 year from the date of purchase as shown by the dealer's stamp and signature and by the purchaser's receipt for tax purposes which must be kept attached to this warranty. The warranty covers only any faults in manufacturing and does not cover accidental damage, incorrect use, negligence or misuse of the product. Use only the accessories provided; the use of accessories other than the approved ones can make your warranty void. Never open the appliance; in the event of any opening or tampering the warranty will automatically become invalid. The warranty does not cover either the parts subject to wear and tear or the batteries(if provided); the warranty expires 1 years from the date of purchase; after this period any technical service shall be changed.

DIAGNOSIS AND TREATMENT

DIAGNOSIS	CHECK POINT	TREATMENT
Unstable value	If you change the place to measure each time	Ref. to "method of correct measurement"
No any display electrodes in screen when press power button.	If wrong electrodes when installing	Please install correct
	If run out of battery	Please charge new battery (LR41)

TECHNICAL SPECIFICATIONS

NAME	Digital thermometer	Measuring range	32°C-42.9°C/90.0-109.9°F
BATTERY	DC1.5V TYPE-LR41	MAX ERROR	0.1°C/0.2°F
RESOLUTION	0.1	WORKING CONDITIONS	10°C-35°C ≤80%RH
STORAGE CONDITIONS	-25°C-55°C ≤85%RH		

CE REP Luxus Lebenswelt GmbH
Kochstr. 1, 47877, Willich, Germany

EC REP Luxus Lebenswelt GmbH
Kochstr. 1, 47877, Willich, Germany

CE 0123

Date of last revision: 24/10/2022 V2

Fabricante/Manufacturer:
Kangfu Medical Equipment Factory
N°380 Ningkang East Road, 325600 Yueqing,
Zhejiang, China

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